



Basic Solution

Alkaline cleaning solution for use in automated analyzers

USES

Basic Solution is employed in automated analyzers as cleaning solution for cuvettes reading, probes and hydraulic system.

PROVIDED REAGENTS

Basic Solution: sodium hydroxide solution, 1 N

INSTRUCTIONS FOR USE

Basic Solution: ready to use.

The solution could be used diluted according to the specifications given for each analyzer.

WARNINGS

The reagent is for "in vitro" diagnostic use.

The solution is irritant. R36/38: irritates eyes and skin. S26/28: in case of contact with the eyes and skin, rinse immediately with plenty of water and seek medical advice. S37/39: wear suitable gloves, Wear eye/face protection.

Use the reagents according to the working procedures for clinical laboratories.

All reagents and samples should be discarded according to the local regulations in force.

STABILITY AND STORAGE INSTRUCTIONS

Basic Solution: stable at room temperature (< 25°C) until the expiration date stated on the box.

INSTABILITY AND DETERIORATION OF REAGENTS

Evident signs of turbidity or presence of particles in suspension indicate reagent deterioration. Discard in such case. Avoid external contamination.

WIENER LAB. PROVIDES









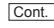


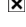

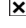
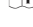





- 1 x 1000 ml (Cat. 1999721).

SYMBOLS

The following symbols are used in packaging for Wiener lab. diagnostic reagents kits.



This product fulfills the requirements of the European Directive 98/79 EC for "in vitro" diagnostic medical devices

-  Authorized representative in the European Community
-  "In vitro" diagnostic medical device
-  Contains sufficient for <n> tests
-  Use by
-  Temperature limitation (store at)
-  Do not freeze
-  Biological risks
-  Volume after reconstitution
-  Contents
-  Batch code
-  Manufactured by:
-  Harmful
-  Corrosive / Caustic
-  Irritant
-  Consult instructions for use
-  Calibrator
-  Control
-  Positive Control
-  Negative Control
-  Catalog number

Wiener Laboratorios S.A.I.C.
Riobamba 2944
2000 - Rosario - Argentina
<http://www.wiener-lab.com.ar>
Dir. Tec.: Viviana E. Cétola
Bioquímica
A.N.M.A.T. Registered product
PM-1102-45



Wiener lab.

2000 Rosario - Argentina